
PROZIS

WHIZZY

Electric Toothbrush

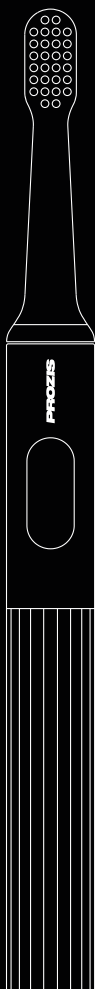User manual
Warranty statement

V1. 12/2019

Prozis is a registered trademark of PROZIS.COM, S.A.

Prozis reserves the right to adjust or modify the product or any of the associated documentation in order to ensure its suitability for use, at any time and without prior notice.

The most recent version of this document, which may contain such adjustments or modifications, is available for reference at www.prozis.com.

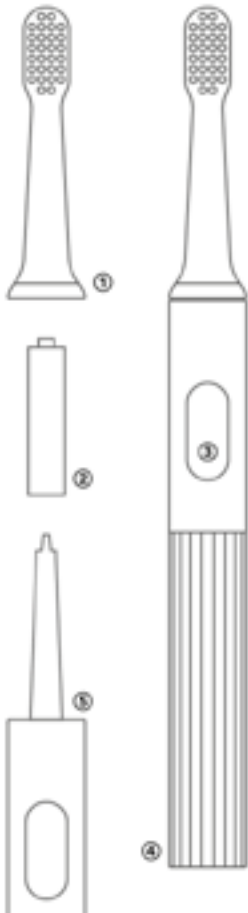
Figure 1

LEGAL COMPLIANCE DECLARATION

By means of this document, PROZIS declares that the **WHIZZY - ELECTRIC TOOTHBRUSH** is fully compliant with all the standards specified in the applicable legislation as well as other relevant clauses.

Prozis will always be available for any additional explanation through our customer support communication channels, found at www.prozis.com/support.

 To avoid any impact on the environment or human health due to dangerous substances present in electrical and electronic devices, the end users of such devices are expected to understand the meaning of the symbol consisting of a crossed-out waste container. Do not dispose of electrical and electronic equipment along with unsorted household waste, instead dispose of it separately and properly.

 Used batteries should not be disposed of along with household waste. They must be placed in a special container for used batteries. This is an obligation intended to ensure the ecological disposal of batteries. Dispose of them only when they are empty.

 The packaging that protects the device against transport damage is made from non-polluting materials that can be disposed of via local recycle bins.

SAFETY WARNINGS

 **CAUTION!** This symbol highlights dangerous situations that can lead to light injuries, or damage, malfunction, and/or destruction of the device.

To prevent harm or injury to users or others, be sure to comply with the below requirements.

- This device can be used by children 8 years old and over, and by people with reduced physical, sensory or mental capabilities or lack of experience and knowledge provided their use is supervised or they have been given instructions concerning the use of the device in a safe way and understand the hazards involved.
- Children must not play with this device.
- Cleaning and maintenance must not be handled by unsupervised children.

LANGUAGE INDEX

EN	3
FR	8
DE	14
IT	20
ES	26
PT	32

INTENDED USE

This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use. Discontinue use of this product and contact a physician/dentist should you experience any discomfort or pain.

MEDICAL WARNINGS

- Please consult your dental professional before you use this product if you have had oral or gum surgery in the last 2 months.
- Please contact your dental professional if excessive bleeding occurs after using this product or should bleeding continue after 1 week of use.
- Please consult your physician before using this product if you have any medical concerns.

WARNINGS

- Do not use the device if it, or parts of it, are visibly damaged.
- Do not alter or modify the device in any way.
- Do not immerse the toothbrush in water or any other liquid.
- Do not operate the device near heat sources.
- Do not use sharp objects to remove dried toothpaste from the handle. Do not puncture the device.

BATTERY WARNINGS

- Do not mix and/or use old batteries and new batteries, or batteries from different brands.
- Do not recharge the batteries.
- Do not throw batteries in a fire.
- Do not connect the + and - electrodes to each other using metal or wire. Do not carry or store batteries along with a metallic objects.
- Avoid direct soldering to batteries.
- Keep batteries away from direct sunlight, high temperatures, and high humidity.
- Check the condition of the batteries on a regular basis. Leaking batteries can damage the device. If the batteries have leaked, use protective gloves. Clean the battery compartment and the batteries with a clean, dry cloth.
- Keep batteries out of reach of small children. Should a child swallow a battery, consult a physician immediately.
- Keep the toothbrush dry before replacing batteries.

CLEANING

- Switch off the device before cleaning it.
- Do not use any chemical or abrasive cleaning agents.
- Do not use sharp objects, nails or brushes.

CLEANING THE TOOTHBRUSH HANDLE

a) Before rinsing the toothbrush body, make sure the battery cover is locked.

b) Remove the brush head and rinse the handle with warm water. Make sure you remove any residual toothpaste. Please do not immerse it in water.

c) Wipe the entire surface of the handle with a damp cloth.

Note: Do not use isopropyl rubbing alcohol, vinegar or bleach to clean handle as this may result in discoloration.

CLEANING THE BRUSH HEAD

a) Rinse the brush head and bristles after each use.

b) Remove the brush head from the handle and

rinse the brush head connection using warm water at least once a week.

REPLACING THE BRUSH HEAD

- Rotate and pull out the used brush head.
- Replace it with the new head and rotate it to lock.

STORAGE

- Store the device away from heat, humidity, and direct sunlight if not in use for a long period of time. To prevent dust and dirt build-up, please use the old packaging to repack the device.
- Remove the batteries if you do not intend to use the toothbrush for a long period of time.

USER MANUAL

PACKAGE CONTENTS

- 1 x Prozis Whizzy - Electric Toothbrush
- 1 x Brush head
- 1 x AAA battery
- 1 x User manual

INSTRUCTIONS

PRODUCT DESCRIPTION *(see page 2)*

1. Brush head
2. AAA battery
3. On/Off button
4. Toothbrush handle
5. Drive shaft

HOW TO USE *(see page 2)*

1. Open the battery compartment by rotating the bottom of the toothbrush handle counterclockwise.
2. Place the batteries according to the polarity specified in the battery compartment.
3. Replace the battery compartment cover, making sure it is securely attached.
4. Rinse the new brush head with fresh water.
5. Press the power button briefly to turn the toothbrush on.

Attaching and Removing the Brush Head

Attaching: Align the brush head with the bristles facing the front of the handle. Match the bars on the brush head with the notches inside the top of the handle and slide the brush head onto the handle.

Removing: Turn the brush head counterclockwise and pull it upwards.

Note: Always use original brush heads to avoid any damage to the toothbrush.

Brushing Technique

- Moisten the brush head under running water.
- Use a pea-sized amount of toothpaste.
- Press the on/off button briefly to switch the toothbrush on.
- Divide your mouth in four parts (quadrants): upper left, upper right, lower left, and lower right. Concentrate on thoroughly cleaning each section as follows:

1. Tilt the brush at a 45-degree angle against the gumline and sweep or roll the brush away from the gumline.

2. Gently brush the outside, inside and chewing surface of each tooth. Move the brush back and forth in circular motions. Clean each quadrant for 30 seconds using the same technique.
3. Gently brush your tongue to remove any bacteria and freshen breath.

- Press the on/off button to switch off the toothbrush.
- Rinse the brush head thoroughly under running water.
- Remove the brush head from the handle. Clean the handle's drive shaft under running water.

TROUBLESHOOTING

Some issues that may arise can be solved by the user. If the problem persists after having tried the following solutions, please contact Prozis customer service.

Technical problem	Cause	Solution
The toothbrush does not work.	The batteries are empty.	Replace the batteries.

WARRANTY STATEMENT

All electrical and electronic products commercialized by Prozis via www.prozis.com are covered by the warranty applicable to the purchase and sale of consumer goods.

GENERAL CONSIDERATIONS

This document contains the terms and conditions of Prozis' warranty for final consumers.

This warranty is valid exclusively for consumers who purchase the product for non-professional purposes.

WARRANTY PERIOD

Prozis guarantees the product is to remain free from material and production defects for a 2-year period, counting from the date of initial purchase by a consumer and respective delivery, or for longer if legally established by the applicable national law.

In the event that repair is required, this 2-year period is suspended during repair time.

CONSUMABLES

The gradual performance decrease caused by prolonged use of the product's consumable components, such as batteries, is not covered by this warranty unless the malfunction is caused by faulty design, materials or production.

HOW TO MAKE A CLAIM

If the need arises to exercise the rights conceded to you through this warranty, please contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions or recommendations. Contact and technical assistance information is also provided at the end of this warranty.

Any claims made regarding this warranty will only be valid if proof of purchase is presented by the original customer. This proof of purchase may be the original invoice or receipt, as long as it in-

cludes information regarding the date of purchase and the product's model name.

WARRANTY REPAIRS

Any repairs made under this warranty cannot be made by third parties. This warranty does not apply to any repairs, or damage caused directly by such repairs, that have impacted or caused the damage that is the subject of the respective warranty claim.

EXCLUSIONS

The following situations are not covered by this warranty:

- Any wear and tear of parts and components resulting from normal product use.
- Cracks, dents, scratches, and other types of superficial damage that only affect the product's appearance.
- Any malfunction caused by inappropriate use of the product that does not comply with user manual instructions.
- Use and storage of the product not complying with user manual instructions.
- Opening of the product by a third party who is not authorized to proceed to its repair.
- Altering or changing parts or components, whether internal or external.
- Damage caused by battery leakage resulting from the breaking of any of its components, or from product misuse.
- Malfunction or damage due to reasons that cannot be ascribed to the production or design of the product.
- Use of the product in the context of a business, occupation or commerce.
- Malfunction due to use that does not comply with the technical or security norms in force, or with user manual instructions.

The services provided by Prozis in order to repair or fix any fault or malfunction resulting from the verification of any of the aforementioned excluded situations will be subject to payment of labor, transport and component costs.

The customer shall bear all risks of loss and damage to the product during transportation to Prozis. This warranty will be void if the returned product arrives with labels or stickers that have been removed, damaged, tampered with or modified in any way.

LIMITED RESPONSIBILITY

Prozis cannot be held responsible for indirect damage or loss due to usage of the product.

WHAT TO DO?

To request warranty service, you must first contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions and recommendations.

Send back the product inside its original packaging, including all accessories and documentation, to the following address:

PROZIS
 Rua do Cais nº 198
 4830-345 Póvoa de Lanhoso
 Braga - Portugal

CONTACTS

Prozis is always available for additional clarification through their customer service communication channel at www.prozis.com/support.

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