
PROZIS

BODY'UP

Electric Muscle Stimulator

User manual
Warranty statement

V.1 11/2019

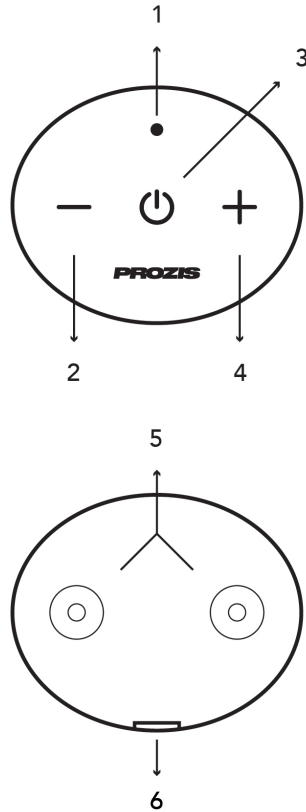Language index	
EN	3
FR	12
DE	22
IT	32
ES	41
PT	50

Prozis is a registered trademark of PROZIS.COM, S.A.

Prozis reserves the right to adjust or modify the product or any of the associated documentation in order to ensure its suitability for use, at any time and without prior notice.

The most recent version of this document, which may contain such adjustments or modifications, is available for reference at www.prozis.com.

Legal compliance statement

By means of this document, Prozis declares that the Body'Up Electric Muscle Stimulator is fully compliant with all the standards specified in the applicable legislation as well as other relevant clauses.

CE Prozis will always be available for any additional explanation through our customer support communication channels, found at www.prozis.com/support.

 To avoid any impact on the environment or human health due to dangerous substances present in electrical and electronic devices, the end users of such devices are expected to understand the meaning of the symbol consisting of a crossed-out waste container. Do not dispose of electrical and electronic equipment along with unsorted household waste, instead dispose of it separately and properly.

 The packaging that protects the device against transport damage is made of non-polluting materials that can be disposed of via local recycle bins.

Safety warnings

 **DANGER!** This symbol, when associated with the note **Danger** means **high risk**. Ignoring this warning may result in injury.

DANGER FOR CHILDREN AND PEOPLE WITH DISABILITIES

- Packing materials are not children's toys.
- Children must not play with the packing material due to risk of suffocation.

- This device is not designed to be used by people (including children) with reduced physical, sensory or mental capabilities or lack of experience and/or understanding, unless supervised by someone responsible for their safety or instructed by them on how to use the device.
- Children must not play with the device.
- Cleaning and maintenance must not be performed by children without supervision.
- Keep the device out of children's reach.

**CAUTION!**

This symbol highlights dangerous situations that can lead to light injuries, or damage, malfunction, and/or destruction of the device.

To prevent harm or damage to users and others, make sure to comply with the following requirements:

- This product is not a medical device. Please refer to the safety warnings before using it.
- If the silicone film in the back of the device is damaged, discontinue its use immediately.
- Do not switch on the device if the product doesn't fit or is not properly in contact with the skin.
- Do not use it for more than 20 consecutive minutes on the same body area.
- Discontinue its use if symptoms such as tingling, itching, redness and inflammation occur.
- If your skin is dry and the device does not fit, you will experience a tingling sensation. You can use a water spray bottle to moisturize your skin.
- Switch the device off after use and before removing it from your body.
- Do not use other electronic products at the same time.

PLEASE DO NOT USE THE DEVICE IN THE FOLLOWING CONDITIONS:

- If you have a cardiac pacemaker, an implanted defibrillator, or any other implant with metal or electronic components.
- If you use an insulin pump.
- If you're running a fever (e.g. >39°C).
- If you have a chronic or acute cardiac arrhythmia, or any disorders of the heart's electrical conduction system.
- If you suffer from any seizure disorders (e.g. epilepsy).
- If you have cancer.
- If you are easily allergic.
- If you have abnormal blood pressure.
- If you have osteoporosis.
- If you're undergoing any treatment.
- If you're pregnant.

**WARNINGS**

- The device must never be used near your heart and must not be placed on any part of the front ribcage, especially on the two large pectoral muscles.
- Do not use across your upper neck and face.
- Do not place over open wounds, sores, or rashes, nor over swollen, red, infected or inflamed areas or skin eruptions (e.g. due to phlebitis, thrombosis or varicose veins).
- Do not place over body areas usually subjected to hair removal, as they are more sensitive.
- Do not use over tattoos.
- Don't use over your abdomen while menstruating.
- Do not place over your carotid sinus nerves, the front of your neck, or around your mouth.
- Do not use in your genital area.
- Keep the device dry. Do not use it in the bathroom or in other damp environments.
- Do not use after drinking alcohol.
- Do not use while driving or operating machinery.
- Do not use while eating.
- Do not use on animals, children or on the speech-impaired.
- Do not use in combination with acupuncture.

**CAUTION!**

- This device is intended for home use only. Do not use it for business purposes.
- Remove your glasses, watches, any jewelry and other metal items from your body before using the device.
- Do not use this device in the presence of electromagnetic fields.
- Do not use this device if the unit has been damaged.
- For hygienic reasons, this device should be used by one person only.
- Switch off the device before attaching or removing it.
- Do not use adhesive tape, bandages or any aid other than the designated gel film to apply the device to your skin.
- Using the device on the same area for a long period of time can irritate the skin. Please consult a dermatologist if the skin irritation persists.
- Pre-treating the skin with greasy creams or lotions is not recommended.

CAUTION! RISK OF MATERIAL DAMAGE

- Protect the product from dirt, moisture and overheating, and use it in dry environments only.
- Do not expose the device to direct sunlight or high temperatures (above 40°C).
- Never submerge the device in water.
- Do not open, disassemble or attempt to repair the device.
- Use this device only with electrodes and accessories recommended by the manufacturer.

MAINTENANCE AND CLEANING

- Switch the device off and remove all plugs before cleaning!
- Do not use chemical or abrasive cleaning agents.
- Do not use scouring milk or steel wool. These will damage the upper surface of the appliance!
- If the main device is dusty, please wipe it using a dry towel.
- Do not wash the device with the accessories.
- Make sure that moisture does not enter the device during cleaning to prevent damage.

HOW TO CLEAN THE ACCESSORIES:

- Separate the device from the accessories.
- Wash separately
- Hand wash only
- Detergents are allowed
- Cold wash (below 30°C)
- Rub using your hands
- Dry in a cool, dry place
- No dry cleaning
- No ironing
- No bleaching
- No soaking
- No machine washing

STORAGE

- Store the device in a cool, dry place, away from direct sunlight.
- Recharge the device on a regular basis if not to be used for a long period of time.
- When the storage period exceeds 6 months, please open the box to check the product and press the power button to see if it's working normally.

User manual

PACKAGE CONTENTS

- 1 × PROZIS Body/Up EMS
- 1 × USB charging cable
- 1 × User manual
- 1 × Spray bottle

PRODUCT DESCRIPTION (SEE PAGE 2)

INSTRUCTIONS FOR USE

Functional area description


- 1 LED light
- 2 Down
- 3 On/Off
- 4 Up
- 5 Magnetic buckle
- 6 Charging port

- The device has 10 adjusted levels of intensity and works for 20 minutes.
- Press the **On/Off** button to switch on the device. The LED light will flash, and the device will beep once.
- Press the **On/Off** button to switch off the device. The LED light will flash, and the device will beep three times.
- After switching on the device, press the **+** or **-** button to adjust the intensity. There are 10 levels of output intensity. When you switch on the device, the default level is the first one. Whenever the intensity is adjusted, the LED light will flash and beep once. When adjusted to the highest or lowest level, the LED light will flash twice in a row as a reminder.
- Should the device be suddenly removed from the skin, the LED light will flash and beep as a reminder. Once reattached, the intensity will be returned to the first level.

Charging Method

- 1 Please use the attached USB cable to charge the device before use. The USB cable can be connected to a computer or adapter and ensures that each charge is enough for 1.5 hours.
- 2 LED light status during charging: the light will flash when the device is charging and will stay on once fully charged.
- 3 The device can be used continuously for about 8 hours after being fully charged.

Accessories

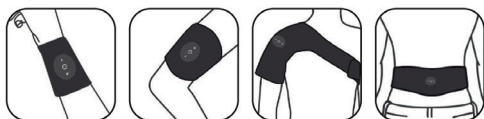
Wrists

Arms

Shoulders

Waist

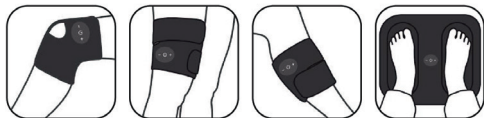
Knees

Thighs

Calves

Soles

Instructions for use

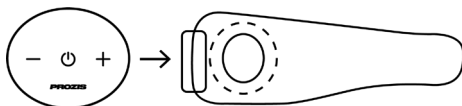- The device and the accessories must be well connected.
- Each accessory has silicone in its back - spraying it with water is recommended.
- Adjust the accessory to your skin.
- Switch it on and the device will start.
- The device works for 20 minutes (automatic program with 3 modes - massaging, tapping and kneading).
- Use the or button to adjust the intensity level (1-10 levels).
- You can switch off the device whenever you feel like it.

Troubleshooting

Some issues that arise may be solved by the user. If the problem persists after having tried the following solutions, please contact Prozis customer service.

PROBLEM	CAUSE	SOLUTION
Stimulation is weak or non-existent	Low intensity	Press the or button to adjust the intensity
	Battery drained	Charge the battery
Device does not switch on	Battery drained	Charge the battery
How to use the device		You should not use it for more than 20 minutes in the same body area.
Where not to use the device		The device should not be used near the heart, over wounds, veins, mouth, genitals, face, and any inflamed or itching areas, or cuts.
Accessory does not fit	Dry skin	If your skin is dry and the device does not fit, you will experience a tingling sensation. You can use a water spray bottle to hydrate your skin.

Warranty statement

All electrical and electronic products commercialized by Prozis via www.prozis.com are covered by the warranty applicable to the purchase and sale of consumer goods.

GENERAL CONSIDERATIONS

This document contains the terms and conditions of Prozis' warranty for final consumers.

This warranty is valid exclusively for consumers who purchase the product for non-professional purposes.

WARRANTY PERIOD

Prozis guarantees the product is to remain free from material and production defects for a 2-year period, counting from the date of initial purchase by a consumer and respective delivery, or for longer if legally established by the applicable national law.

In the event that repair is required, this 2-year period is suspended during repair time.

CONSUMABLES

The gradual performance decrease caused by prolonged use of the product's consumable components, such as batteries, is not covered by this warranty unless the malfunction is caused by faulty design, materials or production.

HOW TO MAKE A CLAIM

If the need arises to exercise the rights conceded to you through this warranty, please contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions or recommendations. Contact and technical assistance information is also provided at the end of this warranty.

Any claims made regarding this warranty will only be valid if proof of purchase is presented by the original customer. This proof of purchase may be the original invoice or receipt, as long as it includes information regarding the date of purchase and the product's model name.

WARRANTY REPAIRS

Any repairs made under this warranty cannot be made by third parties. This warranty does not apply to any repairs, or damage caused directly by such repairs, that have impacted or caused the damage that is the subject of the respective warranty claim.

EXCLUSIONS

The following situations are not covered by this warranty:

- 1 Any wear and tear of parts and components resulting from normal product use.
- 2 Cracks, dents, scratches, and other types of superficial damage that only affect the product's appearance.
- 3 Any malfunction caused by inappropriate use of the product that does not comply with user manual instructions.
- 4 Use and storage of the product not complying with user manual instructions.
- 5 Opening of the product by a third party who is not authorized to proceed to its repair.
- 6 Altering or changing parts or components, whether internal or external.
- 7 Damage caused by battery leakage resulting from the breaking of any of its components, or from product misuse.
- 8 Malfunction or damage due to reasons that cannot be ascribed to the production or design of the product.
- 9 Use of the product in the context of a business, occupation or commerce.
- 10 Malfunction due to use that does not comply with the technical or security norms in force, or with user manual instructions.

The services provided by Prozis in order to repair or fix any fault or malfunction resulting from the verification of any of the aforementioned excluded situations will be subject to payment of labor, transport and component costs.

The customer shall bear all risks of loss and damage to the

product during transportation to Prozis. This warranty will be void if the returned product arrives with labels or stickers that have been removed, damaged, tampered with or modified in any way.

LIMITED RESPONSIBILITY

Prozis cannot be held responsible for indirect damage or loss due to usage of the product.

WHAT TO DO?

To request warranty service, you must first contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions and recommendations.

Send back the product inside its original packaging, including all accessories and documentation, to the following address:

PROZIS

Rua do Cais n° 198

4830-345 Póvoa de Lanhoso

Braga - Portugal

CONTACTS

Prozis is always available for additional clarification through their customer service communication channel at www.prozis.com/support.

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