
PROZIS


PULSE^{hex}

ELECTRIC MUSCLE STIMULATOR

User manual
Warranty statement

V.1 04/2019

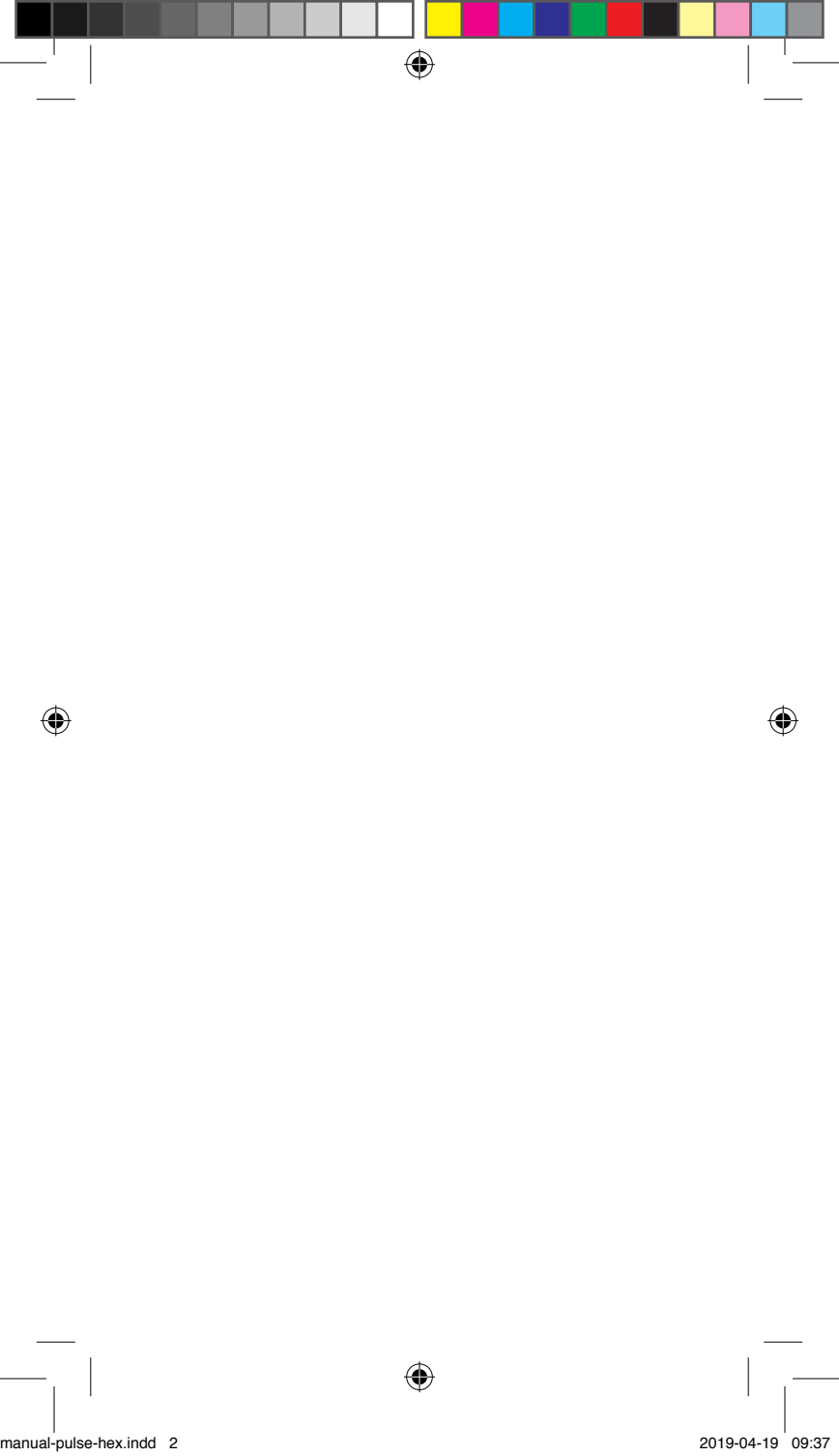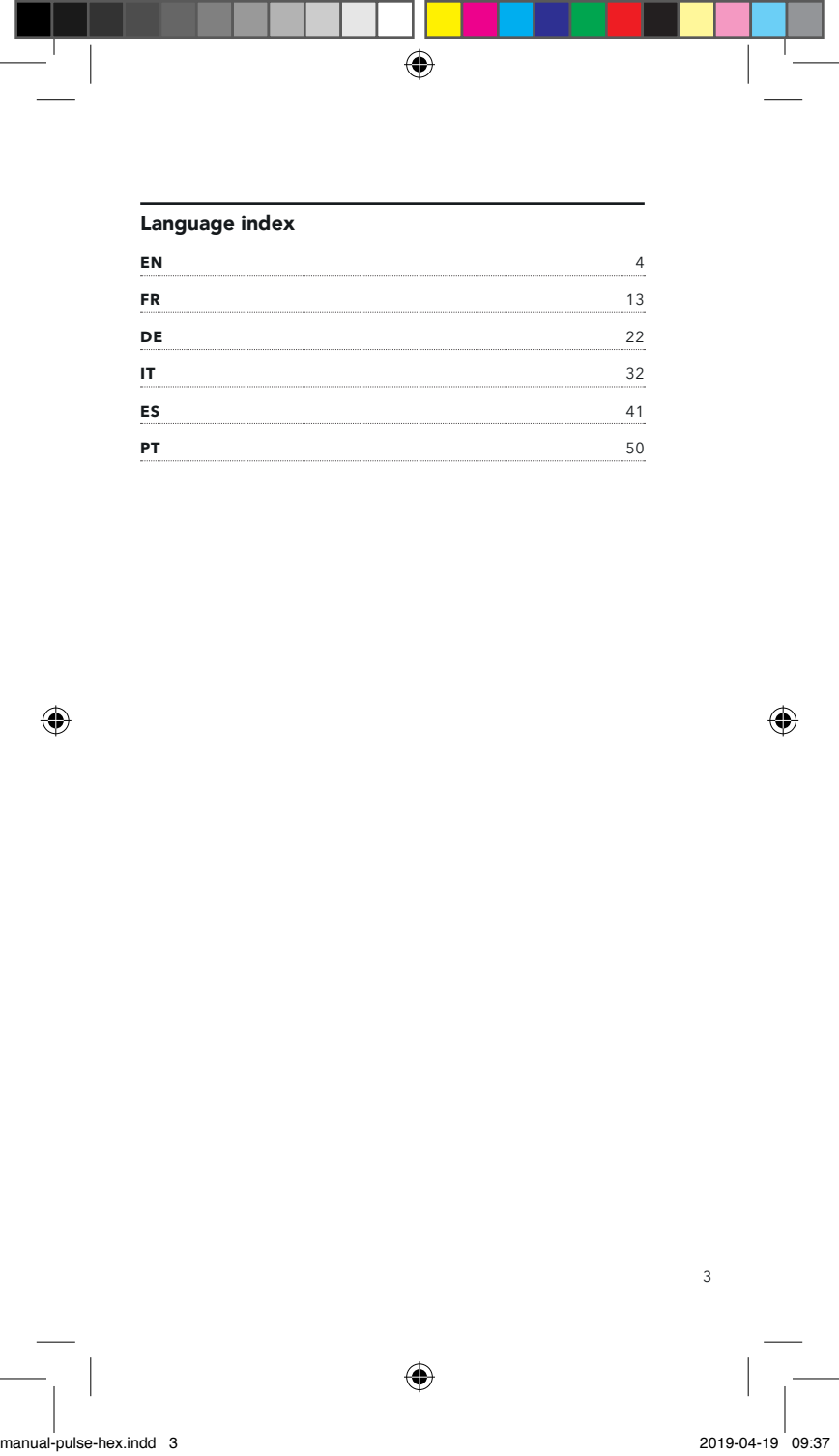Language index

EN	4
FR	13
DE	22
IT	32
ES	41
PT	50

Prozis is a registered trademark of ONSALESIT, S.A.

Prozis reserves the right to adjust or modify the product or any of the associated documentation in order to ensure its suitability for use, at any time and without prior notice. The most recent version of this document, which may contain such adjustments or modifications, is available for reference at www.prozis.com.

Legal compliance statement

By means of this document, Prozis declares that the Pulse-Hex Electric Muscle Stimulator is fully compliant with all the standards specified in the applicable legislation as well as other relevant clauses.

Prozis will always be available for any additional explanation through our customer support communication channels, found at www.prozis.com/support.

To avoid any impact on the environment or human health due to dangerous substances present in electrical and electronic devices, the end users of such devices are expected to understand the meaning of the symbol consisting of a crossed-out waste container. Do not dispose of electrical and electronic equipment along with unsorted household waste, instead dispose of it separately and properly.

® The packaging that protects the device against transport damage is made of non-polluting materials that can be disposed of via local recycle bins.

Safety warnings

CAUTION!

This symbol highlights dangerous situations which can lead to light injuries, or damage, malfunction, and/or destruction of the device.

To prevent harm or injury to users and others, please comply with the following requirements.

- This device is intended for home use only. Do not use for commercial purposes.
- Remove glasses, watches, any jewelry and/or other metal items before using the device.
- Turn the power off after use and before removing the device from your body.
- Do not use in combination with acupuncture.

EN

- Do not use if your hands are wet.
- Keep the device dry. Do not use it in the bathroom or any other humid environments.
- Do not use when charging or replacing the gel pad. Any contact with the device will result in electric shock.
- Keep away from fire and water.
- Never submerge the device.
- Replace the gel pads only when the power is turned off, and with dry hands to prevent device failure, burns or electric shock.
- Do not use across the upper neck and face.
- Keep away from direct sunlight, dust and high temperatures.
- Do not use while driving or operating machinery.
- Do not use while sleeping.
- Do not attempt to disassemble or repair the device yourself.
- Do not use on animals, children or on the speech-impaired.
- The device and battery should be disposed of separately.
- Should you perceive an unusual smell, sound or a response delay during use, turn the device off immediately and contact your reseller.
-  Do not use if you have a cardiac pacemaker, implanted defibrillator, or any other implant with metal or electronic components.
-  The device must never be used near the heart and must not be placed over the front ribcage, especially not on the two large pectorals.
- Do not use in the genital area.
- Do not use during pregnancy.
- Do not place over the carotid sinus nerves, front part of the neck, or over the mouth.
- Do not open the device yourself. It will void the warranty.
- Do not place over or near cancerous lesions.
- Do not place over open wounds, sores, or rashes, nor over swollen,


EN

red, infected or inflamed areas or skin eruptions (due to e.g.: phlebitis, thrombosis or varicose veins).

- Do not use after consuming alcohol.
- Do not use if you have an insulin pump.
- Do not use if you have a fever (e.g. $>39^{\circ}\text{C}$).
- Do not use if you have cardiac arrhythmia, or any heart rhythm or conduction disorder.
- Do not use if you suffer from any seizure disorder (e.g. epilepsy).
- Do not use in the bath or shower.
- Do not use in the presence of electromagnetic fields.
- Do not use this device if the unit has been damaged.
- For hygienic reasons, this device must only be used by one person.
- Do not use adhesive tapes, bandages or any aid other than the designated gel film to apply the device to your skin.
- Stimulation on the same area for a long period of time can irritate the skin. Please consult a dermatologist if the skin irritation persists.
- Pre-treating the skin with greasy creams or lotions is not recommended.

MAINTENANCE AND CLEANING

- Clean the device with a dry or slightly damp cloth. Use only a small amount of soapy water.
- Do not allow the Main Unit to get wet. If this happens, use the device again only once it is fully dried.
- Do not use sharp objects such as nails or hard brush filaments.
- Do not use chemical or abrasive cleaning agents.

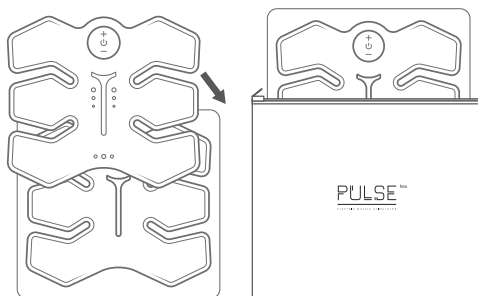STORAGE

- Keep the pad holder clean.
- Store the device in a cool, dry place, away from direct sunlight.
- When not in use, store the device in a safe place to avoid tripping and possible injury.
- Do not store near fire.
- Do not bend the Pad Holder. This will cause deformation that will impede the normal use of the device.


Quick guide & user manual**PACKAGE CONTENTS**

- 1 EMS Control Unit
- 6 Gel Pads
- 1 USB Charging Cable
- 1 Pad Holder
- 1 Carry Bag
- 1 User Manual


EN

INSTRUCTIONS

Functional area description

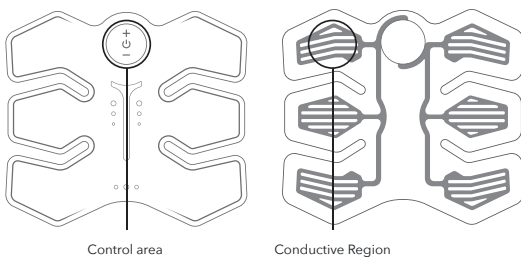Please charge the device using the enclosed USB Cable before use by connecting it to a computer or adapter.

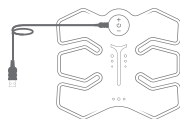- 1.** Remove the blue plastic film from the enclosed electrode Gel Pad.

- 2.** Place the electrode Gel Pad on the back of the control unit.


- 3.** Remove the transparent plastic film and the device is ready to be used.

Please note: This should be done before applying the patch device onto the skin. Avoid wrinkling the stickers. Remove any air bubbles below the sticker surface by pressing them toward the edges.

Buttons

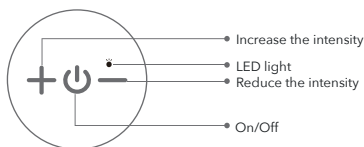- Press the ON/OFF button to turn on the device. It will buzz and the LED light will flash once.
- Press the ON/OFF button again to turn the device off. It will buzz and the LED light will flash three times.
- Press the "+" or "-" button to respectively increase or reduce output intensity. There is a total of 15 intensity levels. After turning it on, the device will operate on intensity level 1. Every time the intensity is adjusted, the device will buzz, and the LED light will flash once. When the highest or lowest intensity is reached, it will buzz twice.

Use

Before using the device, make sure the targeted body parts are clean and dry. At this point, peel the transparent film off the Gel Pads. Place the device onto the targeted area (back, shoulders, arms, legs, etc.) and press down firmly for maximum skin contact. Always start on low intensity – you can gradually increase it according to your preference.

TIPS

- Make sure you replace the battery in time. If you notice weakened buzzing, decreased intensity or a dimmer LED light color, that means the battery is low.
- Turn the power off after use and before removing the device from your body.
- Once removed from the body after use, place the device on the Pad Holder and store it in the Carrying Bag to prevent the gel from getting dirty and drying out.
- The electrode Gel Pads are disposable. When used for 20 minutes each time, they can be used about 30 times, depending on the type and intensity of use. The Gel Pads should be replaced when they have lost their adhesiveness. Keeping the pads clean will extend their lifespan. Avoiding touching the adhesive area with your fingertips will make them last longer.

EN

FREQUENTLY ASKED QUESTIONS

Can I use the device every day?

You can, but do not use it more than once per day on the same body area.

Can it be used continuously?

Do not use it for more than 20 consecutive minutes on the same body area.

What areas of the body should be avoided?

Do not use the device over the heart, veins, face, mouth, genitals, wounds and cuts, moles and birthmarks, or on inflamed or itching skin.

TROUBLESHOOTING

Some of the issues may be solved by the user. If the problem persists after having tried the following solutions, contact Prozis customer service.

TECHNICAL PROBLEM	SOLUTION
The device does not turn on.	Check if the battery is low. It may need to be recharged.
When in use, the device beeps continuously and the intensity cannot be increased.	Make sure the skin is properly clean and the pads are firmly attached to the skin. Also, check if the battery is low and in need of recharging.
The electrode Gel Pads are not adhesive.	Over time, the Gel Pads lose their stickiness. Please contact your reseller for replacements.
The stimulation is too weak or non-existent.	The battery is low and needs to be recharged.

Warranty statement


All electrical and electronic products commercialized by Prozis via www.prozis.com are covered by the warranty applicable to the purchase and sale of consumer goods.

GENERAL CONSIDERATIONS

This document contains the terms and conditions of Prozis' warranty for final consumers. This warranty is valid exclusively for consumers who purchase the product for non-professional purposes.

WARRANTY PERIOD

Prozis guarantees the product is to remain free from material and production defects for a 2-year period, counting from the date of initial purchase by a consumer and respective delivery, or for longer if legally established by the applicable national law. In the event that repair is required, this 2-year period is suspended during repair time.

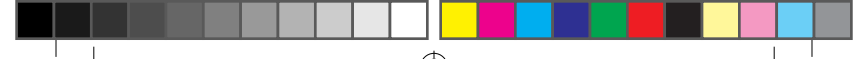EN

CONSUMABLES

The gradual performance decrease caused by prolonged use of the product's consumable components, such as batteries, is not covered by this warranty unless the malfunction is caused by faulty design, materials or production.

HOW TO MAKE A CLAIM

If the need arises to exercise the rights conceded to you through this warranty, please contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions or recommendations. Contact and technical assistance information is also provided at the end of this warranty.

Any claims made regarding this warranty will only be valid if proof of purchase is presented by the original customer. This proof of purchase may be the original invoice or receipt, as long as it includes information regarding the date of purchase and the product's model name.

WARRANTY REPAIRS

Any repairs made under this warranty cannot be made by third parties. This warranty does not apply to any repairs, or damage caused directly by such repairs, that have impacted or caused the damage that is the subject of the respective warranty claim.

EXCLUSIONS

The following situations are not covered by this warranty:

1. Any wear and tear of parts and components resulting from normal product use.
2. Cracks, dents, scratches, and other types of superficial damage that only affect the product's appearance.
3. Any malfunction caused by inappropriate use of the product that does not comply with user manual instructions.
4. Use and storage of the product not complying with user manual instructions.
5. Opening of the product by a third party who is not authorized to proceed to its repair.
6. Altering or changing parts or components, whether internal or external.
7. Damage caused by battery leakage resulting from the breaking of any of its components, or from product misuse.
8. Malfunction or damage due to reasons that cannot be ascribed to the production or design of the product.
9. Use of the product in the context of a business, occupation or commerce.

EN

10. Malfunction due to use that does not comply with the technical or security norms in force, or with user manual instructions.

The services provided by Prozis in order to repair or fix any fault or malfunction resulting from the verification of any of the aforementioned excluded situations will be subject to payment of labor, transport and component costs.

The customer shall bear all risks of loss and damage to the product during transportation to Prozis. This warranty will be void if the returned product arrives with labels or stickers that have been removed, damaged, tampered with or modified in any way.

LIMITED RESPONSIBILITY

Prozis cannot be held responsible for indirect damage or loss due to usage of the product.

WHAT TO DO?

To request warranty service, you must first contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions and recommendations.

Send back the product inside its original packaging, including all accessories and documentation, to the following address:

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CONTACTS

Prozis is always available for additional clarification through their customer service communication channel at www.prozis.com/support.

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