**PROZIS**


**PULSE** hex 2.0

---

ELECTRICAL MUSCLE STIMULATOR

User Manual  
Warranty Statement

V.1 04/2022
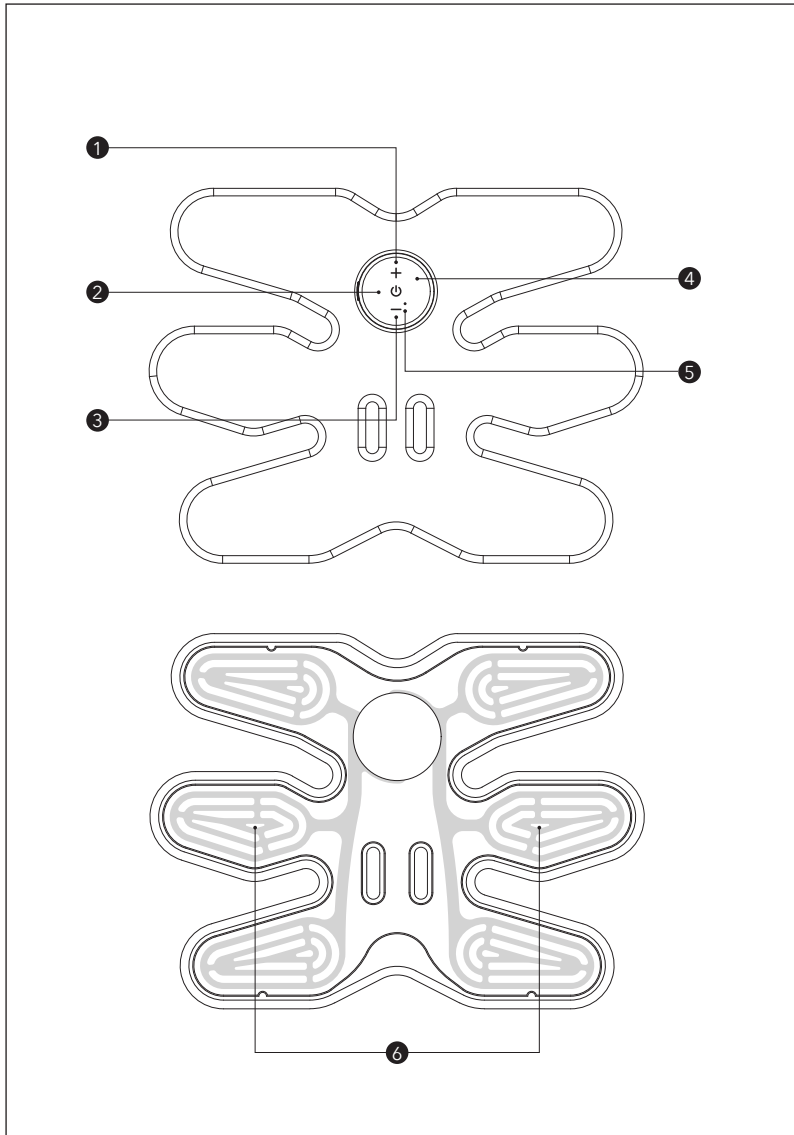
## Language index

<b>EN</b>	4
<b>FR</b>	11
<b>DE</b>	19
<b>IT</b>	27
<b>ES</b>	35
<b>PT</b>	43

Prozis is a registered trademark of PROZIS.COM, S.A.

Prozis reserves the right to adjust or modify the product or any of the associated documentation in order to ensure its suitability for use, at any time and without prior notice.

**NOTE:** The most recent version of this product's user manual can be found at [www.prozis.com/user-manuals](http://www.prozis.com/user-manuals)

## Legal compliance statement

 By means of this document, PROZIS declares that the Prozis Pulse Hex 2.0 - Electrical Muscle Stimulator is fully compliant with all the standards specified in the applicable legislation as well as other relevant clauses.

Prozis will always be available for any additional explanation through our customer support communication channels, available at [www.prozis.com/support](http://www.prozis.com/support)

 To avoid any impact on the environment or human health due to dangerous substances present in electrical and electronic devices, the end users of such devices are expected to understand the meaning of the symbol consisting of a crossed-out waste container. Do not dispose of electrical and electronic equipment along with unsorted household waste. Dispose of it separately and properly instead.

 The packaging that protects the device against transport damage is made from non-polluting materials that can be disposed of via local recycle bins.

## Safety warnings

 **DANGER!** This symbol, in combination with the note "Danger" means a high risk. Failure to observe this warning may result in injury to life and limb.

 **DANGER FOR CHILDREN AND PEOPLE WITH DISABILITIES**

- Packing materials are not toys. Children must not play with the packing material. There is a risk of suffocation.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children without supervision.
- Keep the device out of the reach of children.

 **DANGER! RISK OF ELECTRIC SHOCK!**

- Check the product for damage before use. Do not use the product if damaged.
- Do not modify this device. Do not attempt to disassemble or repair the device yourself.
- Protect your device from impact, moisture, dust, chemicals and sudden temperature

- changes, and keep it away from heat sources (ovens, heaters).
- Keep the device dry. Do not use it in the bathroom or any other humid environments.
- Never submerge the device.
- Do not use it in the bath or shower.

### **DANGER! RISK OF INJURY!**

Do not use if you have a cardiac pacemaker, implanted defibrillator, or any other implant with metal or electronic components.

The device must never be used near the heart and must not be placed over the front ribcage, especially over the two large pectoral muscles.

- Do not use in the genital area.
- Do not use while pregnant.
- Do not place over the carotid sinus nerves, front part of the neck, or over your mouth.
- Do not use across the upper neck and face.
- Do not place over or near cancerous lesions.
- Do not place over open wounds, sores, or rashes, nor over swollen, red, infected or inflamed areas or skin eruptions (e.g., due to phlebitis, thrombosis or varicose veins).
- Do not use if you have an insulin pump.
- Do not use after drinking alcohol.
- Do not use if you're running a fever (e.g., > 39°C).
- Do not use if you have cardiac arrhythmia, or any heart rhythm or conduction disorder.
- Do not use if you suffer from any seizure disorder (e.g. epilepsy).

**CAUTION!** This symbol highlights dangerous situations that can lead to light injuries, or damage, malfunction, and/or destruction of the device.

### **To prevent harm or injury to users and others, please comply with the following requirements.**

- This device is intended for home use only. Do not use it for commercial purposes.
- Remove glasses, watches, any jewelry and/or other metal items before using the device.
- Switch off the device after use and before removing the device from your body.
- Do not use in combination with acupuncture.
- Do not use if your hands are wet.
- Do not use while charging or when replacing the gel pad. Any contact with the device will result in electric shock.
- Keep away from fire and water.
- Never submerge the device.
- Replace the gel pads only when the power is switched off, and with dry hands to prevent device failure, burns or electric shock.
- Keep away from direct sunlight, dust and high temperatures.
- Do not use while driving or operating machinery.
- Do not use while sleeping.
- Do not attempt to disassemble or repair the device yourself.
- Do not use on animals, children or on the speech-impaired.

- The device and battery should be disposed of separately.
- Should you notice an unusual smell, sound or a response while using this device, switch it off immediately and contact your reseller.
- Do not use in the presence of electromagnetic fields.
- Do not use this device if the unit has been damaged.
- For hygienic reasons, this device must only be used by one person.
- Do not use adhesive tapes, bandages or any aid other than the designated gel film to apply the device to your skin.
- Stimulation on the same area for a long period of time can irritate the skin. Please consult a dermatologist if the skin irritation persists.
- Pre-treating the skin with greasy creams or lotions is not recommended.

## User manual

### PACKAGE CONTENTS

- 1 x PROZIS PULSE HEX 2.0
- 1 x Gel Pad Holder
- 6 x Gel Pad
- 1 x USB Charging Cable
- 1 x Carrying Bag
- 1 x Quick Guide

### PRODUCT OVERVIEW (SEE PAGE 2)

1. Increase Intensity
2. ON / OFF
3. Reduce Intensity
4. Control Area
5. LED Light
6. Conductive Region

### MAINTENANCE AND CLEANING

#### Cleaning

- Switch the product off and remove all plugs before cleaning!
- When cleaning the components, do not use aggressive or abrasive cleaning agents or materials such as scouring milk or steel wool. They can damage the upper surface of the appliance!
- Protect the device from water drops and water spray.
- Only use a clean, dry cloth for cleaning it and never use any corrosive liquids.
- Use a slightly damp cloth and a mild cleaning agent to remove any stubborn dirt.
- Do not allow the Main Unit to get wet. If this happens, use the device again only once it is fully dried.
- Make sure moisture does not enter the product while cleaning to prevent damage to the product and related repairs.

#### Storage


- When no longer using the device, store the product in a dry environment, protected from dust and direct sunlight.
- If the device is not going to be used for a long period of time, you should recharge it on a regular basis. This helps maintain battery performance.
- Keep the pad holder clean.
- Do not bend the Pad Holder. This will cause deformation that will interfere with normal use of the device.
- Store your product inside the original case provided.

### INSTRUCTIONS

#### How to use


Before using the device, make sure the targeted body parts are clean and dry.

#### Charging

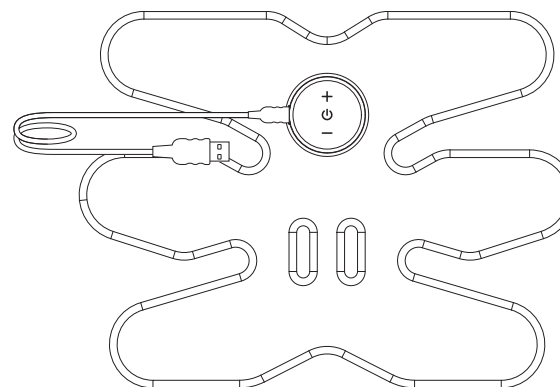Please charge the device using the enclosed USB Cable before use by connecting it to a computer or adapter.


#### Gel Pad Placement

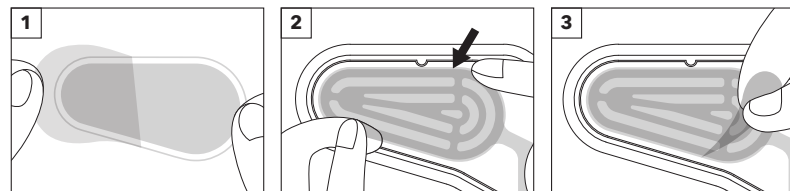1. Remove the blue plastic film from the enclosed electrode Gel Pad.
2. Place the electrode Gel Pad on the back of the control unit.
3. Remove the transparent plastic film and the device is ready to be used.

Please note: This should be done before applying the device onto the skin. Avoid wrinkling the stickers. Remove any air bubbles below the sticker surface by pressing them toward the edges.

### Device Operation

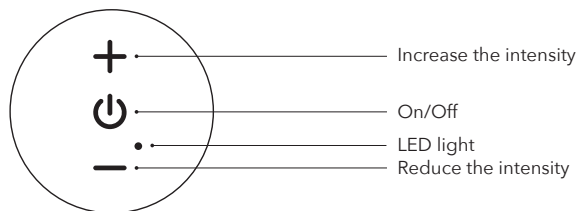- Press the ON/OFF button to switch on the device. It will buzz and the LED light will flash once.
- Press the ON/OFF button again to switch off the device. It will buzz and the LED light will flash three times.
- Press the "+" or "-" buttons to respectively increase or reduce output intensity. There is a total of 15 intensity levels. Once switched on, the device will operate on intensity level 1. The device will buzz and the LED light will flash once every time the intensity is adjusted. It will buzz twice once the highest or lowest intensity has been reached.

Please note: Place the device onto the targeted area (back, shoulders, arms, legs, etc.) and press down firmly for maximum skin contact. Always start on a low intensity - you can gradually increase it according to your preference.

### Tips

- Make sure you charge the battery in time. If you notice weakened buzzing, decreased intensity or a dimmer LED light color, that means the battery is low.
- Switch off the device after use and before removing it from your body.
- Once removed from your body after use, place the device on the Pad Holder and store it in the Carrying Bag to prevent the gel from getting dirty and drying out.
- The electrode Gel Pads are disposable. When used for 20 minutes each time, they can be used about 30 times, depending on the type and intensity of use. The Gel Pads should be replaced when they have lost their adhesiveness. Keeping the pads clean will extend their lifespan. Avoiding touching the adhesive area with your fingertips will make them last longer.

### FREQUENTLY ASKED QUESTIONS

Can I use the device every day?

You can, but do not use it more than once a day on the same body area.

Can it be used continuously?

Do not use it for more than 20 consecutive minutes on the same body area.

What areas of the body should be avoided?

Do not use the device over the heart, veins, face, mouth, genitals, wounds and cuts, moles and birthmarks, or over inflamed or itching skin.

## Troubleshooting

Some issues that arise may be solved by the user. If the problem persists after having tried the following solutions, please contact Prozis customer service.

PROBLEM	CAUSE	SOLUTION
The device is not working.	Battery drained.	Recharge the battery as described.
The stimulation is too weak or non-existent.	Battery drained.	Recharge the battery as described.
When in use, the device beeps continuously and the intensity cannot be increased.	Wrong placement.	Make sure the skin is properly clean and the pads are firmly attached to the skin. Also, check if the battery is low and in need of recharging.
The electrode Gel Pads are not adhesive.	Gel pads lifespan reached.	Over time, the Gel Pads lose their stickiness. Please contact us for replacements.

## Warranty statement

All electrical and electronic products commercialized by Prozis via [www.prozis.com](http://www.prozis.com) are covered by the warranty applicable to the purchase and sale of consumer goods.

### GENERAL CONSIDERATIONS

This document contains the terms and conditions of Prozis' warranty for final consumers.

This warranty is valid exclusively for consumers who purchase the product for non-professional purposes.

### WARRANTY PERIOD

Prozis guarantees the product is to remain free from material and production defects for a 2-year period, counting from the date of initial purchase by a consumer and respective delivery, or for longer if legally established by the applicable national law.

In the event that repair is required, this 2-year period is suspended during repair time.

### CONSUMABLES

The gradual performance decrease caused by prolonged use of the product's consumable components, such as batteries, is not covered by this warranty unless the malfunction is caused by faulty design, materials or production.

## HOW TO MAKE A CLAIM

If the need arises to exercise the rights conceded to you through this warranty, please contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions or recommendations. Contact and technical assistance information is also provided at the end of this warranty.

Any claims made regarding this warranty will only be valid if proof of purchase is presented by the original customer. This proof of purchase may be the original invoice or receipt, as long as it includes information regarding the date of purchase and the product's model name.

## WARRANTY REPAIRS

Any repairs made under this warranty cannot be made by third parties. This warranty does not apply to any repairs, or damage caused directly by such repairs, that have impacted or caused the damage that is the subject of the respective warranty claim.

## EXCLUSIONS

The following situations are not covered by this warranty:

1. Any wear and tear of parts and components resulting from normal product use.
2. Cracks, dents, scratches, and other types of superficial damage that only affect the product's appearance.
3. Any malfunction caused by inappropriate use of the product that does not comply with user manual instructions.
4. Use and storage of the product not complying with user manual instructions.
5. Opening of the product by a third party who is not authorized to proceed to its repair.
6. Altering or changing parts or components, whether internal or external.
7. Damage caused by battery leakage resulting from the breaking of any of its components, or from product misuse.
8. Malfunction or damage due to reasons that cannot be ascribed to the production or design of the product.
9. Use of the product in the context of a business, occupation or commerce.
10. Malfunction due to use that does not comply with the technical or security norms in force, or with user manual instructions.

The services provided by Prozis in order to repair or fix any fault or malfunction resulting from the verification of any of the aforementioned excluded situations will be subject to payment of labor, transport and component costs.

The customer shall bear all risks of loss and damage to the product during transportation to Prozis. This warranty will be void if the returned product arrives with labels or stickers that have been removed, damaged, tampered with or modified in any way.

## LIMITED RESPONSIBILITY

Prozis cannot be held responsible for indirect damage or loss due to usage of the product.

## WHAT TO DO?

To request warranty service, you must first contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions and recommendations.

Send back the product inside its original packaging, including all accessories and documentation, to the following address:

PROZIS  
Rua do Cais n° 198  
4830-345 Póvoa de Lanhoso  
Braga - Portugal

## CONTACTS

Prozis is always available for additional clarification through their customer service communication channel at [www.prozis.com/support](http://www.prozis.com/support)

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