PROZIS

SENSIT^{mini}

smart scale

User Manual
Limited Warranty

EN

Prozis is a registered trademark of PROZIS.COM, S.A.

Prozis reserves the right to adjust or modify the product or any of the associated documentation in order to ensure its suitability for use, at any time and without prior notice.

The most recent version of this document, which may contain such adjustments or modifications, is available for reference at www.prozis.com.

EN

Legal compliance declaration

By means of this document, PROZIS declares that the **SENSIT mini - Smart Scale** is fully compliant with all the standards specified in the applicable legislation as well as other relevant clauses.

Prozis will always be available for any additional clarification through our customer support communication channels, available at www.prozis.com/support.

To avoid any impact on the environment or human health due to dangerous substances present in electrical and electronic devices, the end users of such devices are expected to understand the meaning of the symbol consisting of a crossed-out waste container. Do not dispose of electrical and electronic equipment along with unsorted household waste. Instead, dispose of it separately and properly.

Used batteries should not be disposed of along with household waste. They must be placed in a special container for used batteries. This is an obligation intended to ensure the ecological disposal of batteries. Dispose of them only when they are empty.

The packaging that protects the device against transport damage is made from non-polluting materials that can be disposed of via local recycle bins.

FCC INFORMATION TO THE USER

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.

EN

- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CAUTION!

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

User Manual

The **SENSIT mini** smart scale is a body analysis device that measures weight and uses bioelectrical impedance analysis (BIA) technology to estimate body fat percentage and other body composition data. It is intended for home use only.

WARNING

Please read carefully before use.

- Do not use if you have a pacemaker or artificial lung, or when wearing an ECG or any other internal medical device.
- Do not use while pregnant.
- Do not use on slippery surface, nor when your body is wet.

CAUTION

- Your body composition data measurements may not be accurate when using the scale immediately after intense physical activity.
- This smart scale model is not intended for use by children under the age of 10.
- Smart scales should not be used to diagnose or treat any medical condition. Always consult your physician or healthcare professional to determine what is ideal for you.

TURNING ON THE SCALE

Your smart scale is shipped with 3 AAA batteries. Place the batteries in the battery compartment. Once the scale is powered, it will turn on automatically when you step on it to measure your body weight. However, to measure other body composition parameters such as body fat percentage, the smart scale needs to be used together with a smartphone or tablet with the right software application installed on it.

INSTALLING THE APP AND MEASURING BODY COMPOSITION

The **SENSIT mini** smart scale can connect with smartphones via Bluetooth. Connecting requires the installation of the Prozis Go software application on the respective mobile device. The easiest way to obtain this app is to use your mobile device to scan the following QR code and follow the instructions provided.

EN

EN


Step 1: Getting started

Scan the QR code or search and download the app "Prozis Go" in the App Store or Google Play Store.

Create or use your existing Prozis account.

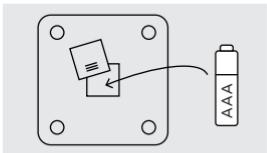
Install the batteries.


Step 2: Pairing the scale with the mobile device

Follow the app installation instructions.

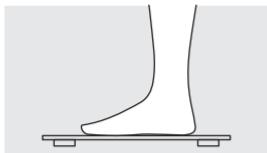
Step on the scale.

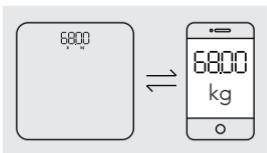Pair the scale with the app.

Step 3: Measuring and checking the analysis report

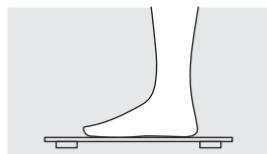Step on the scale again and wait for the number on the display to become stable.

If the number flashes twice, that means the weight measurement is complete. When the number 0 stops moving across the screen from left to right, that means the body analysis is complete.


Check the body analysis report in the app.

MEASUREMENT TIPS


Read the warning on the first page of the manual carefully and follow the instructions. This product employs BIA (bioelectrical impedance analysis) technology. For the most accurate and consistent measurement results, use it during the day and try to maintain the same measuring conditions. **Please avoid measuring:**

- During or immediately after intense exercise
- Immediately after a shower, bath or sauna
- In combination with heavy alcohol consumption
- During or immediately after meals
- While having a fever

When measuring, make sure you stand straight and barefoot on the scale. Wearing socks will lead to inaccurate readings. Dirty feet may also negatively affect measurement accuracy. Via the smartphone app, the measurement units can be changed from imperial/US (pounds for weight, feet/inch for height) to metric (kilograms for weight, meters for height) and vice versa.

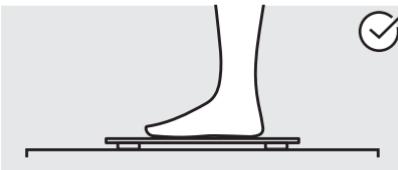Measuring weight correctly: place the scale on a flat, hard surface and reset it before measuring. Before you step on the scale, the weight display should show the value 0.0.

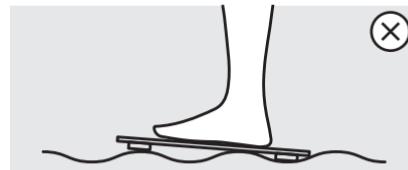Measuring body fat and other body composition parameters correctly: stand straight on the scale, with both feet placed over the electrodes as shown in the above figure.

EN

SPECIFICATIONS


Dimensions:
260 mm × 260 mm × 23 mm
(10.2" × 10.2" × 0.9")

Weight:
1.0 kg (2.2 lb.)

TROUBLESHOOTING

Working normally

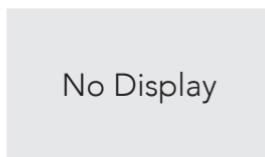Change battery

No Bluetooth connection

Change battery

Overload

Calibrating

WHAT'S INCLUDED

Your **SENSIT mini** smart scale box includes:

- 1 - Smart scale
- 2 - Batteries: 3 x 1.5V AAA
- 3 - User manual

EN

Warranty statement

All electrical and electronic items commercialized by Prozis, via www.prozis.com, are covered by the warranty applicable to the purchase and sale of consumer goods.

GENERAL CONSIDERATIONS

This document contains the terms and conditions of Prozis warranty for final consumers. The warranty is valid exclusively for the consumer who purchase the product for non-professional purposes.

WARRANTY PERIOD

Prozis guarantees that the product is to remain free from material and production defects for a 2-year period, counting from the date of initial purchase by a consumer and the respective delivery, or for longer if legally established by the applicable national law.

In the event that repair is required, this 2-year period is suspended during repair time.

CONSUMABLES

The gradual performance decrease caused by prolonged use of the product's consumable components, such as batteries, is not covered by this warranty unless the malfunction is caused by faulty design, materials or production.

HOW TO MAKE A CLAIM

If the need rises to exercise the rights conceded to you through this warranty, contact the Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions or recommendations. Contact and technical assistance information is also provided at the end of this warranty.

Any claims made regarding this warranty will only be valid if proof of purchase is presented by the original customer. This proof of purchase may be original invoice or receipt, as long as it includes information regarding the date of purchase and the product's model name.

WARRANTY REPAIRS

Any repairs made under this warranty cannot be made by third parties. This warranty does not apply to any repairs, or damage caused directly by such repairs, that have impacted or caused the damage that is the subject of the respective warranty claim.

EXCLUSIONS

The following situations are not covered by this warranty:

EN

1. Any wear and tear of parts and components resulting from normal product use.
2. Cracks, dents, scratches, and other types of superficial damage that only affect the product's appearance.
3. Any malfunction caused by inappropriate use of the product that does not comply with user manual instructions.
4. Use and storage of the product not complying with user manual instructions.
5. Opening of the product by a third party who is not authorized to proceed to its repair.
6. Altering or changing parts or components, whether internal or external.
7. Damage caused by battery leakage resulting from the breaking of any of its components, or from product misuse.
8. Malfunction or damage due to reasons that cannot be ascribed to the production or design of the product.
9. Use of the product in the context of a business, occupation or commerce.
10. Malfunction due to use that does not comply with the technical or security norms in force, or with user manual instructions.

The services provided by Prozis in order to repair or fix any fault or malfunction resulting from the verification of any of the aforementioned excluded situations will be subject to payment of labor, transport and component costs.

The customer shall bear all risks of loss and damage to the product during transportation to Prozis. This warranty will be void if the returned product arrives with labels or stickers that have been removed, damaged, tampered with or modified in any way.

LIMITED RESPONSIBILITY

Prozis cannot be held responsible for indirect damage or loss due to usage of the product.

WHAT TO DO?

To request warranty service, you must first contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions and recommendations.

Send back the product inside its original packaging, including all accessories and documentation, to the following address:

PROZIS
Rua do Cais nº 198
4830-345 Póvoa de Lanhoso
Braga - Portugal

CONTACTS

Prozis is always available for additional clarification through their customer service communication channel at www.prozis.com/support.

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