
**PROZIS**

spot<sup>2</sup>

Vibrating Massager

User manual  
Declaration of warranty  
V.1 03/2020

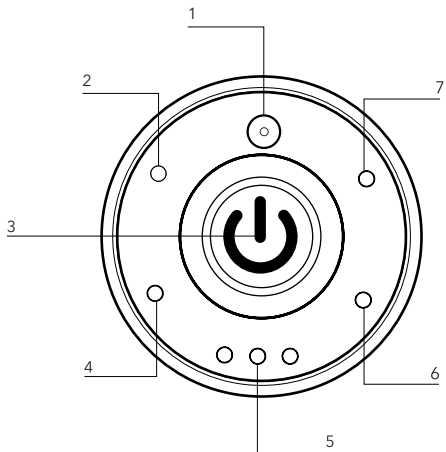---

## Language index

<b>EN</b>	4
<b>FR</b>	22
<b>DE</b>	41
<b>IT</b>	61
<b>ES</b>	79
<b>PT</b>	97

Prozis is a registered trademark of PROZIS.COM, S.A.

Prozis reserves the right to adjust or modify the product or any of the associated documentation in order to ensure its suitability for use, at any time and without prior notice.

The most recent version of this document, which may contain such adjustments or modifications, is available for reference at [www.prozis.com](http://www.prozis.com).

---

## Legal compliance statement

By means of this document, PROZIS declares that the **Prozis Spot<sup>2</sup>** - vibrating massager - is fully compliant with all the standards specified in the applicable legislation as well as other relevant clauses.

 Prozis will always be available for any additional clarification through our customer support communication channels, available at [www.prozis.com/support](http://www.prozis.com/support).

 To avoid any impact on the environment or human health due to dangerous substances present in electrical and electronic devices, the end users of such devices are expected to understand the meaning of the symbol consisting of a crossed-out waste container.

Do not dispose of electrical and electronic equipment along with unsorted household waste. Instead, dispose of it separately and properly.

The packaging that protects the device against transport damage is made from non-polluting materials that can be disposed of via local recycle bins.

---

## Safety warnings

**DANGER!** This symbol, in combination with the note **Danger** highlights a dangerous situation. Failure to observe this warning may lead to limb injury and life risk.

### **DANGER for people with disabilities and children**

- Packing materials are not children's toys. Children must not play with them. There is a risk of suffocation.
- This device can be used by children 8 years old and over, and by people with reduced physical, sensory or mental capabilities or lack of experience and knowledge provided their use is supervised or they

have been given instructions concerning the use of the device in a safe way and understand the hazards involved.

- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children without supervision.
- Keep the device out of the reach of children.

**Caution!** This symbol highlights dangerous situations that can lead to minor to moderate injuries, damage, malfunction, and/or destruction of the device.

**To prevent harm or damage from happening to users or others, make sure to comply with the requirement as below.**

- The device is intended for domestic or private use only, not for commercial use.
- There is a fire risk if the device is used inappropriately or if these instructions are ignored. Therefore, the device must:
  - NEVER be left unattended, especially if there are children nearby.
  - NEVER be covered, such as with blankets, pillows, etc.
  - NEVER be used near petrol or other highly flammable substances.

- The device must not be used close to or underwater (e.g. washbasin, shower, bathtub).
- Ensure that the device is NEVER submerged in water or other liquids. For this reason:
  - NEVER use it outdoors - only in dry indoor room.
  - NEVER place liquid-filled objects, such as pans, near the appliance.
  - NEVER use it in very humid environments.
  - NEVER immerse it in liquids.
- Do not touch the device, the main cable or the main plug with wet hands.
- Before use, make sure the device and accessories display no visible damage. If you have any questions, do not use the device and contact Prozis.
- Never open the housing.
- Keep the device away from pointed or sharp objects.
- Do not use this device if the unit has been damaged.

**This device is intended to massage individual parts of the human body. It cannot replace medical treatment. Do not use the massager if one or several of the following warnings applies to you. If you are unsure whether the massager is suitable for you, please consult your doctor.**

**Do not use this massager:**

- On animals.
- if you suffer from any medical condition or if the area to be massaged is somehow injured (e.g. herniated disc, open wounds).
- For massages in the area around the heart under any circumstances if you have a pacemaker. Should this be the case, please consult your doctor regarding other body parts as well.
- On swollen, burnt, inflamed or damaged skin or body parts.
- Over cuts, capillary or varicose veins, acne, couperose skin, herpes-affected areas or over any other skin diseases.
- While pregnant.
- Over the face (eyes), the larynx or other particularly sensitive body parts.
- While sleeping.
- After taking medication or alcohol (reduced cognition!).
- Inside a vehicle.
- On pubic areas.

**You should consult your doctor before using the massager. This is especially important if you:**

- Have any serious illness.

- Have a pacemaker, an implant or other aids.
- Suffer from thrombosis.
- Have diabetes.
- Are experiencing pain of unknown origin.

The device is only intended for the purpose described in these instructions. The manufacturer is not liable for any damages resulting from inappropriate or careless use.

---

## Cleaning and Maintenance

### CLEANING

- Switch off the device before cleaning!
- When cleaning its components, do not use aggressive or abrasive cleaning agents or materials such as scouring milk or steel wool. These could damage the upper surface of the appliance!
- Protect the device from water drops and water spray.
- Use only a clean, dry cloth to clean the appliance and never any corrosive liquids. Use a slightly damp cloth and a mild cleaning agent to remove any stubborn dirt.

- Make sure moisture does not enter the device while cleaning to prevent damage to the product and any related repairs.
- Do not use the device again until it is completely dry.
- Do not clean the device in a dishwasher.

## **STORAGE**

If you do not plan to use the device for a long period of time, we recommend storing it in the original packaging in a dry environment and ensuring it is not weighed down by additional items.

---

# **User Manual**

## **PACKAGE CONTENTS**

- 1 × Prozis Spot<sup>2</sup> - vibrating massager
- 1 × AC/DC charging adapter
- 1 × User manual

## **SPECIFICATIONS**

Dimensions: 170 (L) x 90 (W) x 90 (H) mm

Weight: 800 g

Voltage: 9 V = 1 A

Power: 9 W

Modes: 3 intensity levels & 1 auto mode

Heating function: up to 30°C

Continuous operating time: 10 min

Battery performance: Up to 5 h\*

Battery: 2200 mAh

Charging time: 2 h

\*Actual battery performance may vary based on settings, environment and usage.

## **PRODUCT DESCRIPTION (SEE PAGE 2)**


1. DC connector
2. Level 3 LED indicator light
3. Multifunction button (Power button)
4. Auto mode LED indicator light
5. Battery LED indicator lights
6. Level 2 LED indicator light
7. Level 1 LED indicator light

**NOTE:** While charging, all three LED lights will flash.

## Instructions

### Body areas over which Prozis Spot<sup>2</sup> can be used


Hips

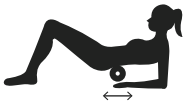Waist

Back

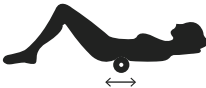Thighs

Arms

Heels

## OPERATION INSTRUCTIONS

### Massage Function

The device provides one auto mode and 3 massage levels, which can be selected by pressing the **Multifunction button**.

#### 1<sup>st</sup> Level

Press the **Multifunction button** and the device will start vibrating at the first level vibration speed. One LED indicator light will switch on. The heating function will also be switched on, and it will take 3-4 min for you to feel the warmth.

#### 2<sup>nd</sup> Level

Press the **Multifunction button** again. Vibration speed will rise to the second level and two LED indicator lights will switch on.

#### Auto Mode

Press the **Multifunction button** for the third time and the device will be set to auto mode. Three LED indicator lights will switch on.

### 3<sup>rd</sup> Level

Press the **Multifunction button** for the fourth time. Vibration speed will rise to level 3 and four LED indicator lights will switch on.

Press the button for the fifth time to switch off the device, as well as the heating function.

On -> Level 1 -> Level 2 -> Auto mode -> Level 3 -> Off

The blue LED indicator lights will switch on according to the selected level.

**Tip:** The device has a built-in 2-second long-press off mode.

This means that once the device is switched on, long-pressing the **Multifunction button** can switch it off, as well as all indicator lights, regardless of its current level.

**NOTE:** Start the device on a low massage level. Discontinue its use should you feel anything out of the ordinary. The device will switch off automatically after 10 minutes. The device should not be used more than 3 times in a row to avoid overheating.

## Battery Charging

- Charge the battery completely before using the device for the first time.
- Plug one end of the supplied charging cable to a regular power socket, and the other end to the device's DC connector.
- Fully charging the battery should take approximately 2 hours.
- The 3 LED indicator lights will flash while charging.
- Once the battery is fully charged, the LED indicator lights will become steady. You may disconnect the power supply.
- Should you feel a drop in massage intensity, stop using the device and recharge it.

**NOTE:** Do not charge the device for more than 6 hours. Unplug the main adapter to avoid danger.

---

## Troubleshooting

Some issues that arise may be solved by the user. If the problem persists after having tried the following solutions, please contact Prozis customer service.

PROBLEM	CAUSE	SOLUTION
The device does not work or stops working while in use.	The charging cable was not properly connected.	Plug in the charging cable correctly.
	The power button did not switch on the device.	Press the power button properly.
	The device has run out of battery.	Charge the battery.
	The device switched off automatically after 10 minutes of continuous use.	Restart the device. Press the power button.
The device is hot or has suddenly stopped working.	You have performed several ten-minute massages.	Remove the power plug and allow the device some time to cool down.

---

## Warranty statement

All electrical and electronic products commercialized by Prozis via [www.prozis.com](http://www.prozis.com) are covered by the warranty applicable to the purchase and sale of consumer goods.

### GENERAL CONSIDERATIONS

This document contains the terms and conditions of Prozis' warranty for final consumers.

This warranty is valid exclusively for consumers who purchase the product for non-professional purposes.

### **WARRANTY PERIOD**

Prozis guarantees the product is to remain free from material and production defects for a 2-year period, counting from the date of initial purchase by a consumer and respective delivery, or for longer if legally established by the applicable national law.

In the event that repair is required, this 2-year period is suspended during repair time.

### **CONSUMABLES**

The gradual performance decrease caused by prolonged use of the product's consumable components, such as batteries, is not covered by this warranty unless the malfunction is caused by faulty design, materials or production.

### **HOW TO MAKE A CLAIM**

If the need arises to exercise the rights conceded to you through this warranty, please contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided

instructions or recommendations. Contact and technical assistance information is also provided at the end of this warranty.

Any claims made regarding this warranty will only be valid if proof of purchase is presented by the original customer. This proof of purchase may be the original invoice or receipt, as long as it includes information regarding the date of purchase and the product's model name.

### **WARRANTY REPAIRS**

Any repairs made under this warranty cannot be made by third parties. This warranty does not apply to any repairs, or damage caused directly by such repairs, that have impacted or caused the damage that is the subject of the respective warranty claim.

### **EXCLUSIONS**

The following situations are not covered by this warranty:

1. Any wear and tear of parts and components resulting from normal product use.
2. Cracks, dents, scratches, and other types of superficial damage that only affect the product's appearance.

3. Any malfunction caused by inappropriate use of the product that does not comply with user manual instructions.
4. Use and storage of the product not complying with user manual instructions.
5. Opening of the product by a third party who is not authorized to proceed to its repair.
6. Altering or changing parts or components, whether internal or external.
7. Damage caused by battery leakage resulting from the breaking of any of its components, or from product misuse.
8. Malfunction or damage due to reasons that cannot be ascribed to the production or design of the product.
9. Use of the product in the context of a business, occupation or commerce.
10. Malfunction due to use that does not comply with the technical or security norms in force, or with user manual instructions.

The services provided by Prozis in order to repair or fix any fault or malfunction resulting from the verification of any of the aforementioned excluded situations will be subject to payment of labor, transport and component costs.

The customer shall bear all risks of loss and damage to the product during transportation to Prozis. This warranty will be void if the returned product arrives with labels or stickers that have been removed, damaged, tampered with or modified in any way.

### **LIMITED RESPONSIBILITY**

Prozis cannot be held responsible for indirect damage or loss due to usage of the product.

### **WHAT TO DO?**

To request warranty service, you must first contact Prozis customer service via the communication channels provided on the Prozis website, and follow the provided instructions and recommendations.

Send back the product inside its original packaging, including all accessories and documentation, to the following address:

PROZIS

Rua do Cais n.º 198

4830-345 Póvoa de Lanhoso

Braga - Portugal

**CONTACTS**

Prozis is always available for additional clarification through their customer service communication channels available at [www.prozis.com/support](http://www.prozis.com/support).

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